

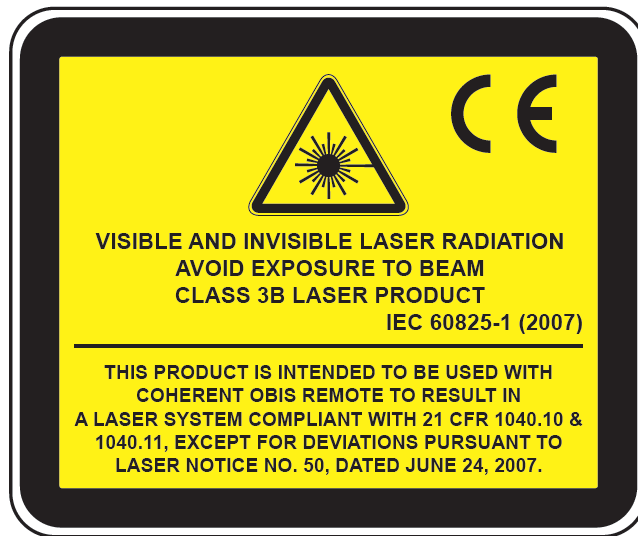


OBIS LX/LS 6-Laser Remote



OBIS LX/LS 6-Laser Remote

The OBIS laser with shutter, OBIS Remote key switch, interlock and 5-second delay allows the OBIS laser system to conform to CDRH Class 3b laser safety requirements.





OBIS LX/LS 6-Laser Remote

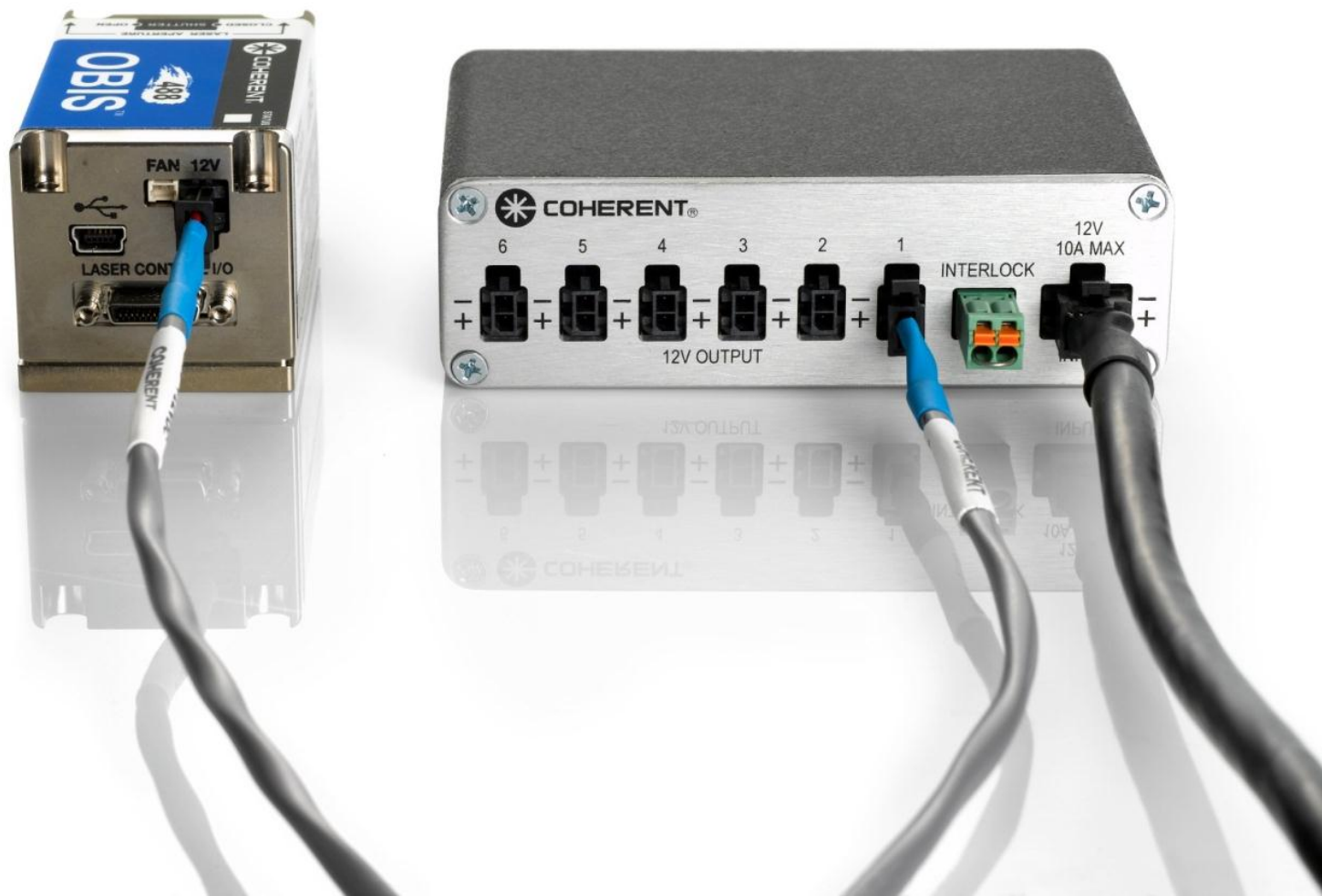
The OBIS 6-Laser Remote can provide an easy installation for power and on/off operation. Includes key switch and interlock.





OBIS LX/LS 6-Laser Remote


Color-coded cables included.





OBIS LX/LS 6-Laser Remote

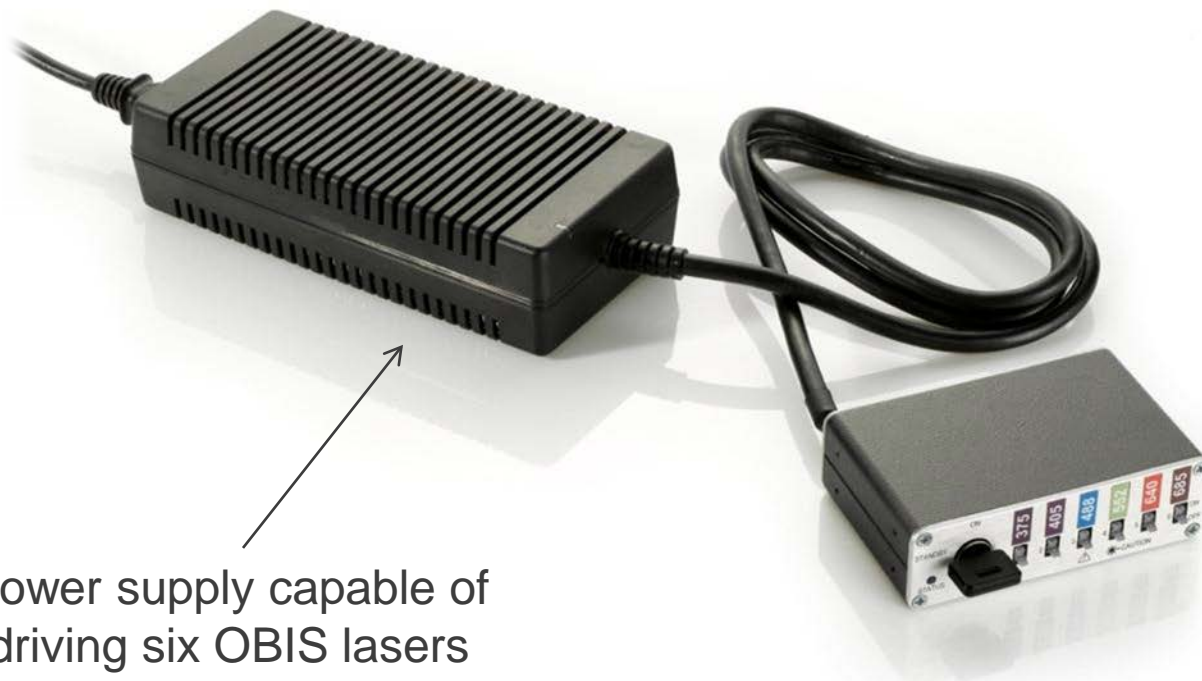Single power supply included to drive six lasers with remote.




Power supply capable of driving six OBIS lasers simultaneously.



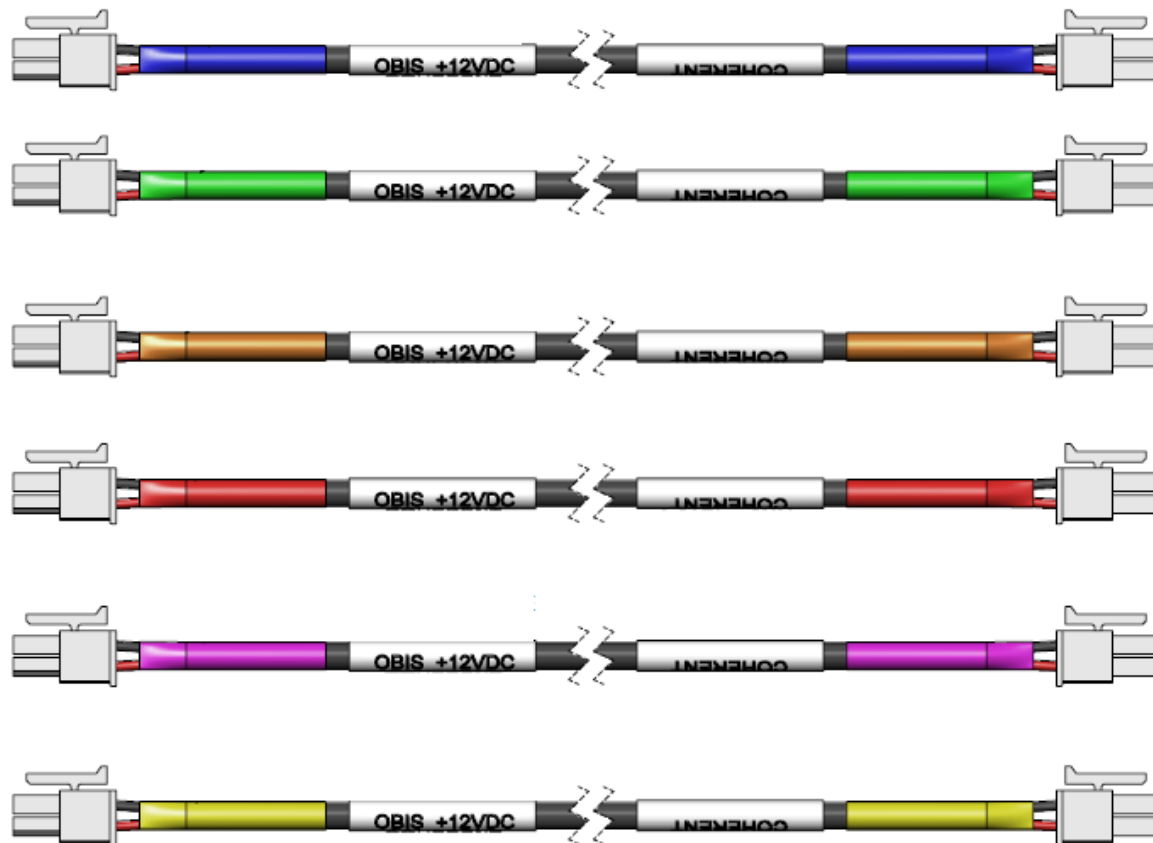OBIS LX/LS 6-Laser Remote


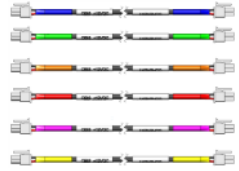

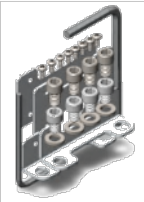





Six remote-to-laser power cables included. Color coded for easy installation and identification. One meter length.





OBIS LX/LS 6-Laser Remote, Part #1203909

The OBIS 6-Laser Remote can add to your laser purchase the following to complete the system.

<p>Laser</p> <p>Not Included</p>	<p>6-Laser Remote</p> 	<p>Laser-to-Remote Cables</p> 	<p>6-Laser Power Supply and Cord</p> 
<p>Bolts and Washers</p> <p>Not Included</p>	<p>Mounting Hardware for Remote</p> 	<p>USB Cable</p> <p>Not Included</p>	<p>Memory Stick includes User Manual</p> 
<p>Laser Safety Pamphlet</p> 	<p>Interlock for Remote</p> 	<p>Keys for Remote</p> 	<p>Wavelength Labels</p> 



OBIS LX/LS Single Laser vs. 6-Laser Remote



OBIS Single Laser Remote

- Full control features and status indicators
- Analog and digital modulation inputs
- USB and RS-232 at the Remote I/O connections



OBIS Six Laser Remote

- Control of up to 6 OBIS Lasers
- Simple power cable to the laser
- NO SDR cable required
- Single Power Supply for 6 lasers